LANAFORM


HEALTH

WF-100


ORAL IRRIGATOR







1

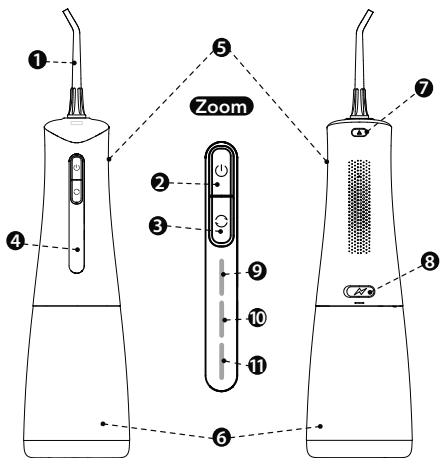2

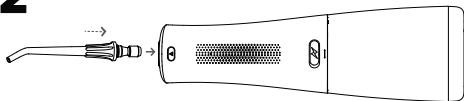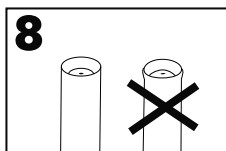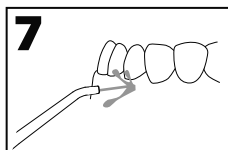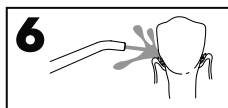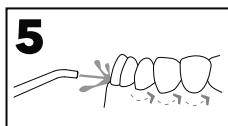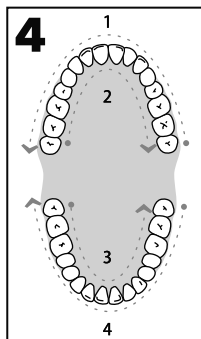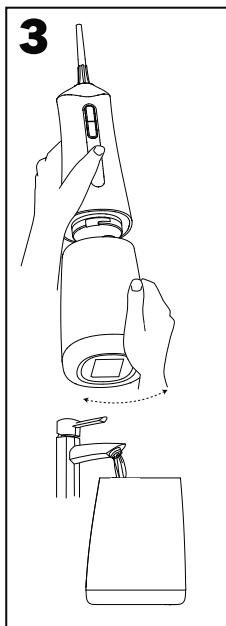ENGLISH

EN

Presentation

Thank you for choosing the LANAFORM oral irrigator.

Thanks to its 4 cleaning modes, the **WF-100** not only cleans the spaces between your teeth by gently dislodging food residues, but also by removing dental plaque and bacteria. In addition, the oral irrigator allows you to gently massage the gums while relieving bleeding.

 Please read all instructions before using your product, particularly the basic safety requirements.

 The photographs and other product images in these instructions and on the packaging are as accurate as possible but do not claim to be a perfect reproduction of the product.

Contents

- 01 Unpacking
- 02 Package contents
- 03 Symbols
- 04 Safety precautions
- 05 Product description
- 06 Instructions for use
- 07 Cleaning and maintenance
- 08 Changing the tip
- 09 Technical specifications
- 10 Troubleshooting
- 11 Warranty
- 12 Disposal of waste

1. Unpacking

Please open the package carefully before use. Check if all accessories are available and if any components have been damaged during transport. Follow the installation and operating instructions in these instructions for use. In the event of any damage or malfunction, please contact the reseller or LANAFORM directly. In the event of a complaint, you will need to have the device model, serial number, date of purchase, your contact details and your address.

2. Package contents

- 1 oral irrigator
- 4 tips
- 1 USB charging cable
- User manual

3. Symbols

The following marks and symbols may appear in the instructions for use.

Symbol	Description
	This symbol warns you of the risk of injury or health hazards
	Prohibited action.
	Fire hazard
	Note
	Consult instructions for use.
	CE sign This product meets the requirements of the applicable European and national directives
	Electrical waste must be sent to a dedicated collection point for recycling.
	Use only indoors

4. Safety precautions

Read the following precautions carefully before using the oral irrigator.

Caution

- **Before starting to use the device, read the safety instructions carefully.** Failure to follow the following instructions may result in physical injury or material damage. Keep this user manual and make it available to all other users. If you pass on the device to anyone else, give them these instructions for use as well.
- **Do not direct the water jet towards the back of the tongue, the inside of the ear, nostrils or other delicate areas.**
- **This device is capable of producing pressures that could cause serious injury in these areas.**
- **This device is for personal use only and is not intended for commercial use.**
- **This device is not intended to be used by people, including children, with reduced physical, sensory or mental capabilities, or people with lack of experience or knowledge, unless they are supervised or instructed on the use of the device by a person responsible for their safety.** Children should be supervised to make sure they do

not play with the device.

- **Oral irrigators can be used under supervision by children at least 8 years of age.**
- **The device is designed only for the use described in these instructions for use.** Use the device only as directed. The manufacturer cannot be held responsible for damage caused by inappropriate or improper use.
- **Use only the tips and accessories recommended by LANAFORM.**
- **Do not drop or insert any objects into the holes of the device.**
- **Do not use or operate the product outside, in a place where aerosol cans are used or in a place where oxygen is administered.**
- **Do not use iodine, saline solution (salt solution), essential oils, mouthwash or any other chemical solution in this product's tank.** The use of these products may reduce the performance of the device and shorten its service life.
- **Remove any oral jewellery items before using this product.**
- **Please ensure that the appearance of the product is not**

damaged and that the start light is normal before using this product.

- **People with periodontal disease, who have undergone dental treatment, or who are concerned about symptoms inside their mouth should consult a healthcare professional before use.**

- **If, during use or after using the device, you notice any injury to your gums, teeth or mouth, stop using the device and consult your healthcare professional.**

- **Do NOT use the device if you have an injury to your mouth, gums or teeth.**

- **Oral irrigators should be used with cold or tepid water (40°C maximum).**

- **Use the device only for oral hygiene. This could cause an accident or injury if the water jet is directed towards the eyes, nose or throat.**

- **Do not share the tip with your family or others. This may cause infection or inflammation.**

- **Do not dry the main unit with a hairdryer or a fan heater. This may cause malfunction or**

deformation of the parts.

- **Do not wipe with diluents, benzene, alcohol etc. This could cause malfunction, cracking or discolouration of the parts.**

- **Store the device in a location with a low humidity level after use. Leaving the device in a bathroom may cause malfunction.**

Electrical shock

- **As with any electrical device, the oral irrigator should be used with care and caution to avoid any electric shock.**

- **Make sure to use a charger with a DC 5V 1A output that complies with national standards when charging this product and properly connect the charger plug to the power outlet to avoid any malfunction or fire caused by heat.**

- **Only plug the adapter into the socket when the device is turned off.**

- **Never use the device if it or its accessories are visibly damaged.**

- **Do not use the device if it has fallen into water.**

- **Do not use the device during**

a storm.

- If the device is faulty or breaks down, switch it off immediately. Disconnect the plug from the socket.

- Do not pull on the power cable or on the device to disconnect the mains adapter from the socket.

- Never hold or carry the device by its power cable.

- Do not pinch the cable.

- Do not use the oral irrigator without first installing the tip.

- Since this product has an integrated rechargeable battery, do not place the product in a fire and avoid using it in a high temperature and overheated environment, as a malfunction or fire caused by heat may occur.

- Keep the mains adapter away from hot surfaces.

- Make sure that the power cable and mains adapter plug do not come into contact with water, steam or other liquids.

- Disconnect the power plug from the domestic socket when not in use for an extended period of time. Otherwise, this could cause an electric shock

or fire due to electrical leakage resulting from a deterioration in insulation.

- Never pick up a device that has fallen into water. Immediately disconnect the mains power adapter from the socket.

- Unplug the device before filling the tank.

Repair

- The device must not be opened.

- Only a specialist can repair electrical appliances. Improper repair may result in a significant hazard to the user. For any repairs, contact customer services or an authorised dealer. Do not attempt to service the oral irrigator unless you are a qualified professional. Only qualified maintenance professionals are authorised to carry out internal maintenance if necessary.

Fire hazard

- Any improper use or use that does not comply with the instructions for use results in a fire hazard!

- Do not use the device in a covered area, under a blanket or

cushions, for example.

- **Do not use the device near flammable materials.**
- **Make sure that the air vents are not covered.**
- **Never disconnect the device when switched on. Always turn it off first.**
- **Do not expose the device to high temperatures or an excessively humid environment.**

5. Product description

- 01 Tip
- 02 On/Off switch
- 03 Mode selection button
- 04 Battery / Operating mode indicator LED
- 05 Main unit
- 06 Water tank
- 07 Tip fixation/ release switch
- 08 USB cable connector
- 09 High mode
- 10 Medium mode
- 11 Soft mode

6. Instructions for use

Charging

- 01 Please charge the product before using it for the first time. Make

sure to use the charger according to point 4. If you use it abroad, first check whether the country's voltage and charger are suitable for this product.

- 02 Switch off the product first when charging. Connect the charger to the USB cable that is included. Insert the charger into the mains socket and connect the product to the USB charging cable.
- 03 The battery indicator shows the charging status. Once charging is complete, the indicator will remain static. Press and hold the mode selection button to display the battery status. At this point, disconnect the product from the USB cable and disconnect the charger from the power supply. It's quite normal for the product and charger to be slightly warm.
- 04 Four hours are more or less needed to fully charge the product using the charger that meets the requirements.

The time may vary depending on the specific circumstances.

- 05 If the battery indicator flashes while the device is in use, it tells the user to charge the device. When the battery is empty, the indicator will flash 4 times and then automatically turn off.
- 06 The product cannot be used while it is being charged.
- 07 To extend the battery life, charge the device immediately when the battery level is low.

Warning

Do not insert or remove the power cord with wet hands (this may cause electric shock and short circuit).

Do not use it in places where water is accessible, such as in a bath or in a humid place, to avoid malfunction or short circuit.

You must use the charger that meets the charging requirements.

You should not be too close to the product during charging; when you leave for a long

period, you must stop charging the device.

Please immediately turn off power when faults occur during charging. If the USB cable is damaged or the plug is loose, please stop charging immediately.

Please use and store the USB cable correctly. Do not fold, stretch, twist, press or wrap around the charger.

How to use

Firmly attach the tip to the instrument using the fixation switch. 

Before filling the tank, please make sure that the device is switched off. Open the water tank by turning the main unit counter-clockwise, keep the tank upright and fill the tank with water ensuring that the temperature of the tank does not exceed 40°C. 

Screw the water tank back onto the main unit in the clockwise direction and make sure it is securely locked.

To turn on the device, hold the on/off switch for 1 second to turn on the device.

Choose one of the cleaning

modes by successively pressing the mode selection button.

There are 4 cleaning modes:

- "Soft" mode

This mode is intended for new users and people with sensitive teeth. To apply this mode, select the bottom light on your device.

- "Medium" mode

This mode is intended for daily cleaning of dental residues, suitable for most users. Apply this mode by selecting the light in the middle of your device.

- "High" mode

The powerful pressure in this mode is intended to clean up more stubborn dental residues and dental plaque. Select the upper light on your device to apply this mode.

- "Pulse" mode

Alternating the pressure of the sprayed water, this mode provides a gentle massage of the gums while removing dental residues. To select this mode, press the mode switch several times until all lights turn on alternately.

Once your treatment has been carried out, press the on/off button 2 times to switch off the

device. For safety reasons, the device has an automatic shut-down after 2 minutes. If you want to use it again, check the water level of the tank and fill it if necessary before pressing the on/off switch.

Directions for use: 4

For best results, start with the molars and move the device forward in the mouth to finish with the incisors. During cleaning, position the tip at a 90° angle to the gum and move it along the gingival line. Stay between each tooth for a short time. Close your lips slightly to prevent water from splashing out, but leave enough space for the water to flow freely from your mouth to the sink.

Interdental cleaning 5

Direct the water jet to the space between the teeth. This action allows you to thoroughly clean the interdental space that is inaccessible with a toothbrush.

Cleaning the periodontal pockets 6

Direct the water jet into the space between the teeth and the gum and slowly move the water jet along the gum.

Use soft mode if your gums are sensitive.

Gum care 7

Stimulate your gums with the massage provided by the water jet. Use soft mode if your gums are sensitive.

For each of these steps, carefully pass the water jet over the front and back of the teeth and behind the molars. This will ensure optimal cleaning.

Device memory

The device will remember the last mode used for your dental treatment. This mode will be activated directly the next time you use it.

After each use, open and empty the tank. Replace it and drain the remaining water by turning on the device until the oral irrigator no longer projects water at all. Then switch it off (as indicated above in the instructions for use).

Note

Place the tip in your mouth and gently close your mouth before pressing the on/off switch.

For use, place your face above a sink and, with the device

upright and tip installed, lift your elbow.

Pressing the mode selection button during use will change the mode.

Turn the tip to adjust the direction of the water flow.

The first time you use it, start with "soft" mode. The water pressure can be gradually increased so that the oral cavity adapts. You may also notice a little discomfort in your gums during the first few uses. This is normal and should stop after 1-2 weeks. However, if this persists or if you experience abnormal pain, stop using the oral irrigator immediately and consult a healthcare professional.

After use, press the on/off switch to turn the device off.

Note

- 01** The device might not operate at an ambient temperature of around 5°C or less.
- 02** Do not push the cannula tip too hard against the teeth or gums.
- 03** Since the plaque adheres to the surface of the teeth cannot be removed only with water, brushing is always necessary.

It is recommended to clean the inside of your mouth with water after brushing your teeth.

- 04** Do not run the main unit on empty except to empty the remaining water after use. Using the main device without water may cause a malfunction.
- 05** If the tip of the cannula is deformed, replace it with a new one recommended by LANAFORM.  **8**

7. Cleaning and maintenance

Cleaning the main unit:

- Only clean the device with clean water or a neutral detergent. Do not use corrosive or abrasive detergents. Do not use water above 50°C on any part of the device as this could cause damage or loss of performance.
- Do not immerse the main unit in water. If necessary, clean the surface with a soft cloth and non-corrosive detergent.

Tank cleaning:

- Detach the tank from the main unit, as explained in point 6.

- Rinse the tank with lukewarm water. Do not clean the reservoir with detergent of any kind.
- Replace the tank after cleaning.

Limescale cleaning:

- Depending on the mineral content of the water used, a limescale deposit may appear on the device. This may reduce the performance of the device.
- Point the tip towards the sink, turn on the device and leave it on until there is no water in the tank. Replace the tank with clean hot water (maximum 40°C) and repeat the above-mentioned tank cleaning steps. This cleaning process should be repeated every 1-3 months in order to maintain product performance.

8. Changing the tip

- The tip is a consumable. The tip should be replaced approximately every 6 months for hygiene reasons, even if its shape remains unchanged.  **8**

9. Technical specifications

Product Name	WF-100
Product reference	LA140304
Dimensions	22.4 x 7 x 7.9 cm
Weight	+/- 285 gr
Battery specification	2000 mAh
Charging time	+/- 4 hrs
Power consumption excluding charge	Maximum 5W
Charger requirement	5V, 1A
Standby current	Less than 20uA

10. Troubleshooting

 Problem	 Cause	 Solution
There is a leak between the tip and the main unit	The tip holder seal is damaged or the tip has not been correctly installed	Contact the seller for repair or reinstall the tip correctly.
There is a leak from the water injector	The seal is damaged or there are impurities in the injector	Contact the seller for repair

 Problem	 Cause	 Solution
There is a leak in the water tank	Unit closure system is not watertight	Replace seal or contact seller for repair
Water pressure is insufficient	The device is not charged enough	Try again after the device is fully charged
The device does not start	The device is not charged enough	Try again after full charge of the device or contact the seller for repair
The device does not spray water	The device is used for the first time or has not been used for a certain period of time	Add enough water to the tip installation hole and try again after reinstalling the tip

11. Warranty

EN LANAFORM guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The LANAFORM warranty does not cover damage caused by normal product wear and tear. In addition, the warranty for this LANAFORM product does not cover any damage caused as a result of misuse or improper use, accidents, attaching any unauthorised accessories, modification of the product or any other condition of any kind beyond the control of LANAFORM.

LANAFORM shall not be held liable for any type of incidental, consequential or special damage.

All implied warranties of fitness for the product are limited to a period of two years from the original date of purchase, provided that a copy of the proof of purchase can be presented.

Upon receipt, LANAFORM will repair or replace your device, as appropriate, and return it to you. A warranty claim can only be made through the LANAFORM Service Centre. Any maintenance of this product entrusted to any person other than the LANAFORM Service Centre cancels this warranty.

12. Disposal of waste

The packaging consists entirely of environmentally safe materials that can be deposited at your local sorting centre to be used as secondary materials. Cardboard can be placed in a paper collection container. Packaging films should be dropped off at your local sorting and recycling centre.

If you are no longer using the device, dispose of it in an environmentally friendly manner and in accordance with legal guidelines.

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