


Recommended Cameras for the “DELTA 20” using HEINE SLR Photo Adaptor

Manufacturer	Type	HEINE SLR Photo Adaptor
Canon	EOS 7D / 50D / 60D	K-000.34.185
	EOS 350D*/400D*	
	EOS 450D/500D/550D	
	EOS 600D/1000D/1100D	
Nikon	D 40* / D 80* / D 40X	K-000.34.186
	D 60 / D 90	
	D 200 / D 300S	
	D 700 / D 7000	
	D 3000 / D 3100	
	D 5000 / D 5100	
Olympus	E-450 / E-510	K-000.34.187
	E-520 / E-620	

* no longer available

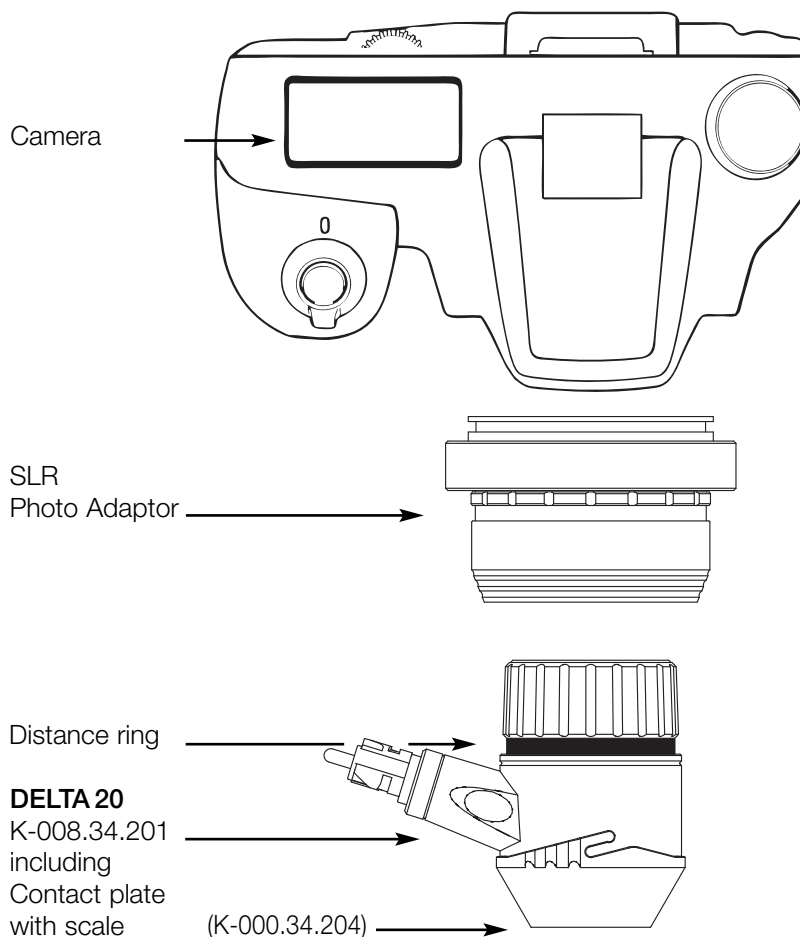

Camera setup for SLR-Adaptor

Nikon

- ISO: 400
- Program: M
- Exposure: 1/40s
- White balance: manual

Canon

- ISO: 400
- Program: Tv
- Exposure: 1/40s
- White balance: manual

Olympus

- ISO: 400
- Program: S
- Exposure: 1/40s
- White balance: manual

Canon Camera Name Chart

Europe	USA	Japan	MPixel
Canon EOS 300D	Digital Rebel	Kiss Digital	6
350D	Rebel XT	Kiss N	8
1000D	Rebel XS	Kiss F	10
400D	Rebel XTi	Kiss X	10
450D	Rebel XSi	Kiss X2	12
500D	Rebel T1i	Kiss X3	15
550D	Rebel T2i	Kiss X4	18